

Achilles™ NC

0.018" system

Combination of strength and optimal control in routine and challenging cases

Broad range of sizes available

Confident to treat all lower extremity lesions

Non-Compliant OTW

Excellent for high pressure dilatation and lesion preparation

Seamless tip-to-shaft design


Provides smooth transmission and maximized pushability

Low tip profile 0.020"

Facilitates lesion entry

Optimized non-compliant balloon material

High pressure balloon inflation up to 22atm

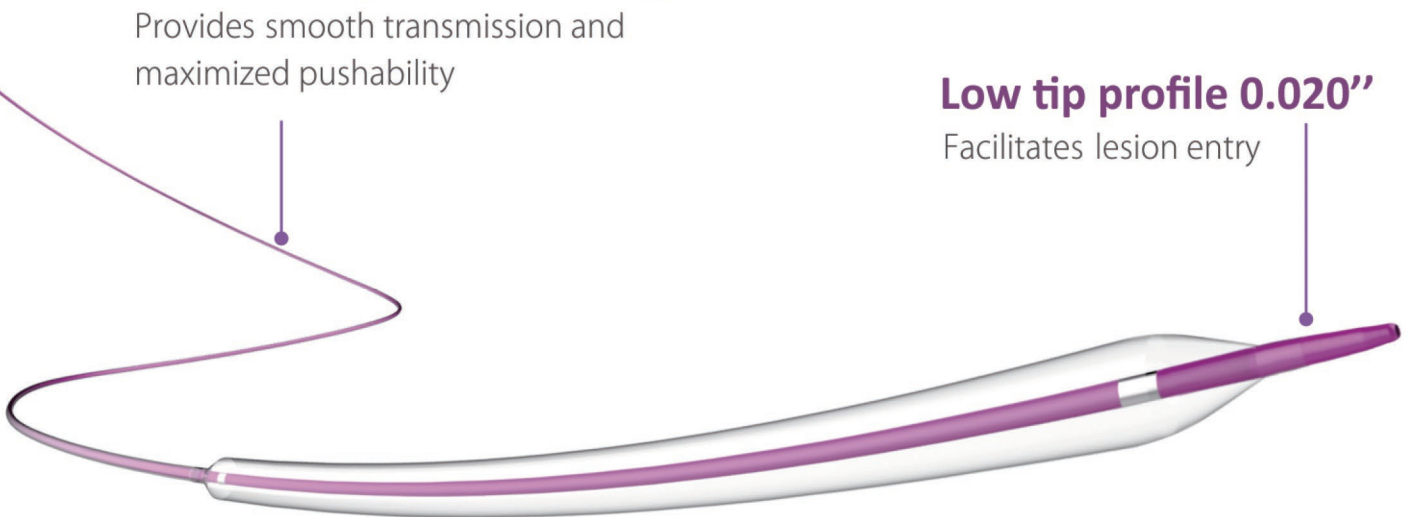

Specifications

Balloon Diameter (mm)	Balloon Length (mm)									Catheter Working Length (cm)	Sheath Compatibility (F)	Nominal Pressure (atm)	Rated Burst Pressure (atm)
	10	15	20	30	40	60	80	100	150				
3.0	—	—	V	V	V	V	—	—	—	40	4	12	22
3.5	—	—	V	V	V	V	—	—	—	40	4	12	22
4.0	—	—	V	V	V	V	—	—	—	40	4	12	22
4.5	—	—	V	V	V	V	—	—	—	40	5	12	20
5.0	—	—	V	V	V	V	—	—	—	40	5	12	20
5.5	—	—	V	V	V	V	—	—	—	40	5	12	20
6.0	—	—	V	V	V	V	—	—	—	40	6	12	20
7.0	—	—	V	V	V	V	—	—	—	40	6	12	18
8.0	—	—	V	V	V	V	—	—	—	40	6	12	18
2.0	V	V	V	V	V	V	V	V	V	70	4	12	22
2.5	V	V	V	V	V	V	V	V	V	70	4	12	22
3.0	V	V	V	V	V	V	V	V	V	70	4	12	22
3.5	V	V	V	V	V	V	V	V	V	70	4	12	22
4.0	V	V	V	V	V	V	V	V	V	70	4	12	22
4.5	V	V	V	V	V	V	V	V	V	70	5	12	20
5.0	V	V	V	V	V	V	V	V	V	70	5	12	20
5.5	V	V	V	V	V	V	V	V	V	70	5	12	20
6.0	V	V	V	V	V	V	V	V	V	70	6	12	20
7.0	V	V	V	V	V	V	V	V	V	70	6	12	18
8.0	V	V	V	V	V	V	V	V	V	70	6	12	18
2.0	V	V	V	V	V	V	V	V	V	90	4	12	22
2.5	V	V	V	V	V	V	V	V	V	90	4	12	22
3.0	V	V	V	V	V	V	V	V	V	90	4	12	22
3.5	V	V	V	V	V	V	V	V	V	90	4	12	22
4.0	V	V	V	V	V	V	V	V	V	90	4	12	22
4.5	V	V	V	V	V	V	V	V	V	90	5	12	20
5.0	V	V	V	V	V	V	V	V	V	90	5	12	20
5.5	V	V	V	V	V	V	V	V	V	90	5	12	20
6.0	V	V	V	V	V	V	V	V	V	90	6	12	20
7.0	V	V	V	V	V	V	V	V	V	90	6	12	18
8.0	V	V	V	V	V	V	V	V	V	90	6	12	18
2.0	V	V	V	V	V	V	V	V	V	150	4	12	22
2.5	V	V	V	V	V	V	V	V	V	150	4	12	22
3.0	V	V	V	V	V	V	V	V	V	150	4	12	22
3.5	V	V	V	V	V	V	V	V	V	150	4	12	22
4.0	V	V	V	V	V	V	V	V	V	150	4	12	22
4.5	V	V	V	V	V	V	V	V	V	150	5	12	20
5.0	V	V	V	V	V	V	V	V	V	150	5	12	20
5.5	V	V	V	V	V	V	V	V	V	150	5	12	20
6.0	V	V	V	V	V	V	V	V	V	150	6	12	20
7.0	V	V	V	V	V	V	V	V	V	150	6	12	18
8.0	V	V	V	V	V	V	V	V	V	150	6	12	18



BrosMed Medical Co., Ltd.

Add: 15th building, SMEs Park, Songshan Lake Hi-Tech Zone, 523808, Dongguan, Guangdong, China
 Tel: +86 (769) 2289 2018 Fax: +86 (769) 2289 2016
 E-mail: sales@brosmed.com www.brosmed.com

BrosMed Medical B.V. (International Sales Office)

Mgr. Buckxstraat 8, 6134 AP Sittard The Netherlands
 Tel: +31 (0) 850 140 901
 E-mail: cs@brosmed.com
 Web: www.brosmed.com

US FDA K160941



GRA-6021/Rev.01